


SCANX[®] 14 PORTABLE DIGITAL IMAGING SYSTEM

INTRODUCING A FAST AND AFFORDABLE
DIGITAL SYSTEM YOU CAN TAKE WITH YOU.

IF YOU ARE A PRACTITIONER WHO REQUIRES RADIOGRAPHS ON THE GO, THE SCANX 14 PORTABLE OFFERS A NEW WORLD OF EASE AND CONVENIENCE. IMAGINE THE TIME AND EFFICIENCY YOU WILL ADD TO EACH DAY WITH A PORTABLE DIGITAL IMAGING SYSTEM THAT ENABLES YOU TO VIEW IMAGES ON YOUR LAPTOP IN SECONDS FOR ON-THE-SPOT DIAGNOSES.

- ✓ ScanX 14 Portable rapidly produces diagnostic quality images that are ready for viewing in just 55 seconds on your laptop. Acquire images, accurate the first time. No need for multiple exposures, "just in case."
- ✓ Use your existing x-ray equipment to expose a reusable Phosphor Storage Plate, the same way you expose film.
- ✓ Never again deal with film, chemistry, processor maintenance or duplicating. No need for a darkroom. ScanX 14 operates in normal room light.
- ✓ Enhance images with supplied compatible imaging software to maximize the diagnostic information.
- ✓ ScanX technology saves you time and space with easy digital file storage and archiving. Retrieve and review patient x-rays, send them to a colleague or burn a copy on CD at the push of a button.
- ✓ ScanX 14 Portable facilitates easy mobility and is equipped with a power inverter for use in your vehicle.

ScanX 14 is sold, installed and serviced through ALLPRO IMAGING's long standing nationwide network of authorized dealers. To locate your authorized ALLPRO IMAGING dealer, call direct: **800.247.8324 ext 5611.**

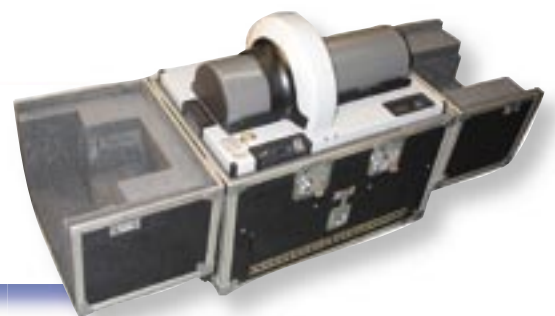

 **ALLPRO IMAGING**
THE IMAGING EXPERTS

SCANX 14 PORTABLE

SPECIFICATIONS

Dimensions (open case):	61" W x 20" D x 27" H
Dimensions (storage):	31" W x 20" D x 28.5" H
Weight:	113 lbs. (in case)
Electrical:	115V, 60 Hz
Imaging Media:	Phosphor Storage Plates 14"x17", 10"x12", 8"x10"
Laser Classification:	Class 1 Laser Product Compliance with FDA HHS 21 CFR 1040.10 and IEC 60825-1
System Includes:	Case, Eraser, Power Supply & Computer Connector Cord, three sets of plates, plate protectors and cassettes
Software:	Included

COMPUTER SYSTEM REQUIREMENTS/ RECOMMENDATIONS

CPU Speed:	2.8 GHz Pentium 4 Windows 2000 Service Pack 4, or XP Pro Service Pack 1a or later
RAM:	1 GB
Hard Drive:	120 GB
CD ROM Drive:	48X
Monitor:	2 MP 20.1" Class Monochrome LCD or SVGA Monitor (1024x786 resolution)
Video:	32 MB RAM
Keyboard, Mouse:	Standard
Backup Device:	Recommended for daily saving
External Surge Protector:	Recommended
USB Port:	Must be 1.1 or higher
Printer:	Optional

ACCESSORIES



PSP's (shown with Plate Protectors) are available in 8" x 10", 10" x 12" and 14" x 17" sizes.



Cassettes. Available in 8"x10", 10"x12" and 14"x17" sizes.



The included ScanX Eraser thoroughly and uniformly erases PSPs.

JUST RIGHT FOR THE OFFICE

SCANX 14 (with built-in Eraser)

Load plates and read images in just 55 seconds. The ScanX 14 has a built-in eraser. The eraser uniformly and thoroughly erases the Phosphor Storage Plates.

SPECIFICATIONS

Dimensions:	22" D x 24" W x 64" H	System Includes:	Power Supply and Computer Connector Cord
Weight:	120 lbs.	Software:	Included
Electrical:	115V, 60 Hz	Computer Requirements:	See ScanX 14 Portable
Imaging Media:	Phosphor Storage Plates 14"x17", 10"x12", 8"x10"		



 **ALLPRO IMAGING**
THE IMAGING EXPERTS

www.allproimaging.com

Contact your authorized ALLPRO IMAGING dealer or call direct: 800.247.8324 ext 5611
70 Cantiague Rock Road, P.O. Box 870, Hicksville, New York 11802 • 516.433.7676 • Fax 516.433.7683